

# SHOWER TROLLEYS



Also Featuring | Change Tables | Mortuary Tables

Manufactured in Australia by

 **MODSEL**  
careful by design

Distributed in South Australia by

  
**SURGICAL AND MEDICAL SUPPLIES PTY LTD**

## introducing...

Modsel Shower Trolleys are constructed for both patient safety and ease of use by allied health professionals, and feature a clean design with durable materials that are ideal for harsh environments. They are manufactured in Australia so support, service, repairs and spare parts are readily available.

These shower trolleys can be used to shower or bath patients including bariatric patients.

The trolley design allows for easy patient transfers with fold away dropsides and fold-down liner walls.

Liners are made of heavy duty, non-stick vinyl and are replaceable and repairable. The liner can be filled with water when the drain hose is hung up on the end rail.

They include a pillow, drain outlet and drain hose.

The trolley foot end can be lowered to improve drainage.



### **Aquaenius**

Standard model  
180Kg capacity  
Hydraulic operation or fixed height  
Available as a 'Junior' size

Specifications : page 4



### **Aquatuff**

Bariatric model  
300Kg capacity  
Hydraulic or electric operation  
Available as a 'Junior' size  
Wide or extra wide top and liner options available

Specifications: see page 6

### **Aquatuff** electric backrest, split dropsides



### **Aquatuff** extra wide



### **Aquatuff** Junior top



## popular options and applications

Both styles of shower trolley can be supplied with the liner replaced by a mattress for use as a Change Table.

“Tailor the purpose of your trolley”



Both styles of shower trolley can be supplied without a liner for use as a Mortuary Trolley.

A concealment cover is also available as an option.



# ► aquarius

standard shower trolley

## special features

- 1 Economical hydraulic trolley
- 2 Concealment cover option Mortuary Trolley version
- 3 Junior size top option

“Modsel design for budget constraints.”



## summary

We are able to bring our Modsel shower trolley design quality to you at an affordable price, helping to meet tight budget constraints and maximising the number of working units available.

The refined framework and proven mechanics awards long life value and consists of considered, easy to clean surfaces coupled with large castors to ensure the highest level of manoeuvrability.

Our shorter top option on this model creates flexibility to cater for smaller sized patients and space limitations.

## standard features

- Foot end can be lowered for improved drainage
- Liner wall folds down to cover any gap between bed and trolley
- Liner repairable with Clean Patch Vinyl
- Can fill with water to use as a bath when hose hung on end rail
- Includes pillow

## technical specs

### Aquarius shower Trolley: ST- AQM

Length - 1950mm (+120mm length adjustment)  
Junior Length - 1680mm (+120mm length adjustment)  
Width -790mm  
Liner Base - 1800x560mm  
Castor Size -75mm  
3 x Break Wheels + 1 Steer Wheel  
Height Range - 680 to 880mm hydraulic  
Fixed Height version available  
Max Load 180kg  
Construction - Triple coated steel with Stainless Steel side and end rails  
Standard Liner Colour - Royal Blue


MADE IN AUSTRALIA

## custom liner colours



## popular options

- 1 Custom liner colours
- 2 Change table mattress
- 3 Bariatric Aquatuff with electric back rest

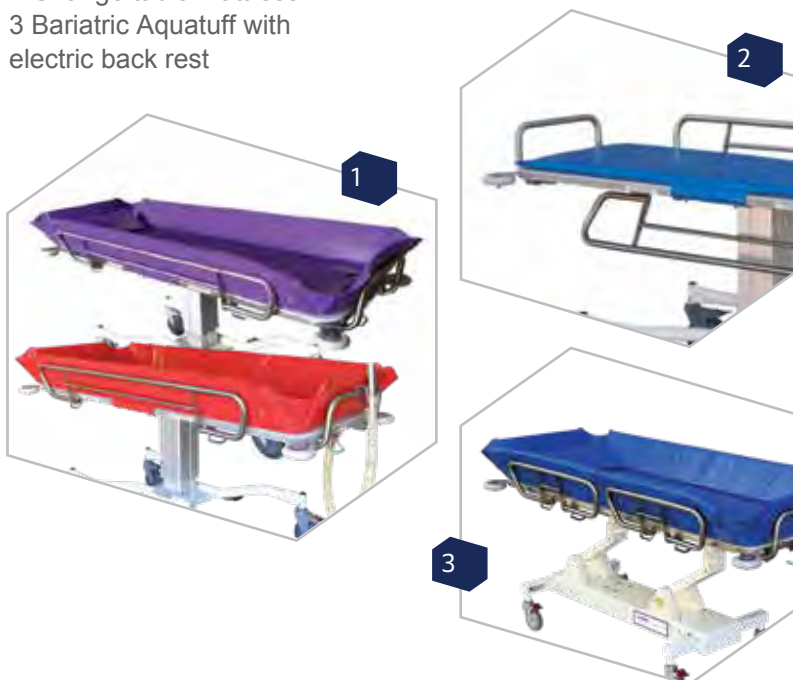

## more options

Custom hose lengths  
Change Table (mattress replaces shower liner)  
Mortuary Trolley (no mattress or shower liner)



We take the harsh wet room environment seriously. The Modsel Shower Trolley range is designed to withstand cleaning products and moisture over a long period of time – every part is considered to assist infection control, from the double coated inner frame to fully welded, outer stainless steel dropsides and end rails.

Cleaning is quick and easy with the removable bath liner, constructed from heavy duty non-stick vinyl which protects the patient's skin whilst enduring for many years the demanding task of every day bathing and showering.

Transfer patients safely with the zero gap dropsides and twin release catches preventing accidental operation. All patient types are accommodated for with shorter, longer and wider top options, bariatric and electric height adjustment with powered backrest available.

# ► aquatuff

bariatric shower trolley

## special features

- 1 Electric or hydraulic height adjustment
- 2 Electric backrest option
- 3 Bariatric wide top options

“Robust and completely reliable.”



## summary

Having equipment that is robust and completely reliable is of utmost importance when caring for larger or heavier patients.

We have increased componentry strength in the design and the decision to use total stainless steel reinforces the lifetime value that this shower trolley will give.

The hydraulic or manual choice operates to safely lift 300kg without overload, whilst four twin wheels with brakes, wider top and electric backrest options all increase versatility in special needs.

## standard features

- Foot end can be lowered for improved drainage
- Electric Backrest includes split dropsides
- Liner wall folds down to cover any gap between bed and trolley
- Liner repairable with Clean Patch Vinyl
- Can fill with water to use as a bath when hose hung on end rail
- Includes pillow

## technical specs

### Aquatuff Shower Trolley: ST- AQT

Length - 1950mm (+120mm extra length adjustment)  
Width - 790mm  
Liner Base - 1800 x 560mm  
Bariatric Width - 100or 400mm wider  
Castor Size - 125mm  
3 x Brake Wheels + 1 Dual Brake/Steer Wheel  
Height Range - 620 to 900mm  
Max Load - 300kg  
Construction - Fully Stainless Steel  
Standard Liner Colour - Royal Blue

MADE IN AUSTRALIA

## custom liner colours



## popular options

- 1 Custom liner colours
- 2 Change table mattress
- 3 Economical Aquarius trolley



## more options

Custom hose lengths  
Change Table (mattress replaces shower liner)  
Mortuary Trolley (no mattress or shower liner)



We take the harsh wet room environment seriously. The Modsel Shower Trolley range is designed to withstand cleaning products and moisture over a long period of time - every part is considered to assist infection control, from the double coated inner frame to fully welded, outer stainless steel dropsides and end rails.

Cleaning is quick and easy with the removable bath liner, constructed from heavy duty non-stick vinyl which protects the patient's skin whilst enduring for many years the demanding task of every day bathing and showering.

Transfer patients safely with the zero gap dropsides and twin release catches preventing accidental operation. All patient types are accommodated for with shorter, longer and wider top options, bariatric and electric height adjustment with powered backrest available.

Contact | Chris Lee | 0476 777 767



**SURGICAL AND MEDICAL SUPPLIES PTY LTD**

16 Kensington Road Rose Park SA 5067

Telephone: (08) 8332 1666 Fax: (08) 8332 5132

**Email: [sales@surmed.com.au](mailto:sales@surmed.com.au)**

**Website: [www.surmed.com.au](http://www.surmed.com.au)**

Disclaimer: The information contained on this brochure is for general information purposes only. Surgical and Medical Pty Ltd accepts no responsibility for errors or omissions in the contents. For further clarification on any aspect of the information contained herein, please refer to the technical data sheets from the manufacturer.