

SA distributor




SURGICAL AND MEDICAL SUPPLIES PTY LTD
16 Kensington Road Rose Park SA 5067
Telephone: (08) 8332 1666 Fax: (08) 8332 5132
Email: sales@surmed.com.au
Website: www.surmed.com.au

ACTION[®]

BODY SMART[®] POSITIONING

Pressure Care Products

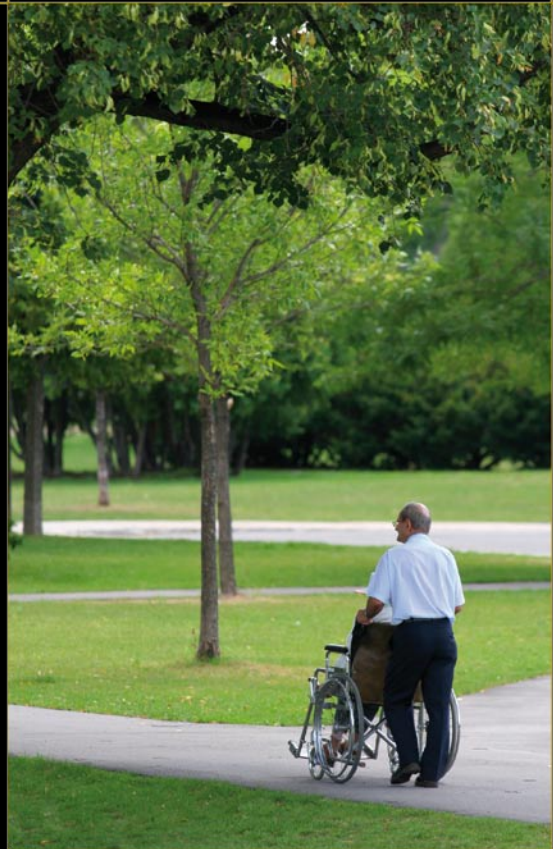

***For pressure protection
....everywhere***

Action® Mattress Overlays


Item #	Size	Weight
6100	Mattress Overlay 69 x 43 x 2.2 cm	6.5 k
6300	Mattress Overlay 68.6 x 116.9 x 2.2 cm	17.5 k
6301	Full Mattress Overlay 68.6 x 203.2 x 1.58 cm	22 k
6301H	Full Mattress Overlay w/Hitch Same as 6301	22 k
6302	Full Mattress Overlay 88.9 x 203.2 x 1.58 cm	30 k
6302H	Full Mattress Overlay w/Hitch Same as 6302	30 k
6303	Full Mattress Overlay 68.6 x 188 x 1.58 cm	20.3 k
6303H	Full Mattress Overlay w/Hitch Same as 6303	20.3 k
6304	Shear Reduction Overlay 68.6 x 188 x .8 cm	12.2 k
6400	Combo Overlay 40.6 x 85 x 2.2 cm (unfolded) 40.6 x 40.6 x 4.4 cm (folded in cover)	8 k
6400H	Combo Overlay w/Hitch Same as 6400	8 k
6405	Combo Overlay 40.6 x 95.2 x 2.2 cm (unfolded) 45.7 x 40.6 x 4.4 cm (folded in cover)	9 k
6405H	Combo Overlay w/Hitch Same as 6405	9 k



Made of 100% Akton® polymer, these Mattress Overlays are used to cover the entire surface or isolated areas of an existing standard mattress.

Used by hospitals and long term care facilities as a significant component in prevention and treatment of decubitus ulcers, AKTON polymer will not migrate like gels, does not promote bacterial growth, and is ideal for multiple client use. Some models are available with a hook and loop hitch used to secure one end of the overlay to the head of the bed.

The hitch will reduce slippage when the bed is raised.

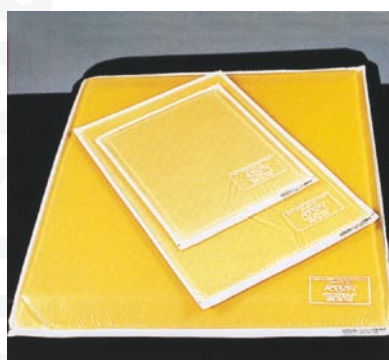

Action® Small Pad

Model 40402 27.9 x 22.9 x 0.95 cm

Action® Medium Pad

Model 40403 43 x 27.9 x 0.6 cm

This pad is larger than the model 40402. Suitable uses include: padding for the iliac crest, under heels and elbows while patient is in lateral position. Fits well within a pillow slip to relieve pressure on the head.

Action® Large Pad

Model 40404 56 x 50 x 1.3 cm

Action® Commode Pads

Action® Commode seat pads provide comfort and relief from pressure and shear. Hook and loop make product easy to attach, use, and maintain.

Item#	Size
COM1616	40.6 x 40.6 x 1.59 cm Oval Opening
COM1818	45.7 x 45.7 x 2.54 cm Split Opening



Action[®] Heel Protector Pad

Model 40401

32 x 20 x 1.3 cm

This heel protector pad should be used for any procedure that is lengthy and involves pressure at the patient's heel area. It is easy to use. Care providers simply place patient's foot in the pad and secure hook and loop attachment.




Action[®] Heel Support

Model 40502

15 x 8.3 x 4.5 cm

Provides the same protection as model 40401, however, configured to secure and protect heel area without hook and loop attachment. Place patient's heel in cup cut-out. Slit in back of cup secures the patient's Achilles tendon area.

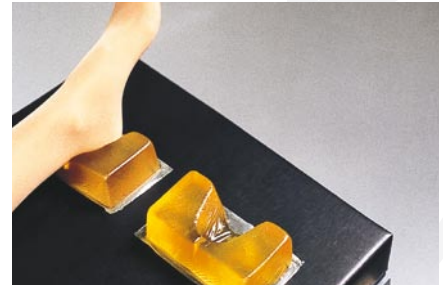

Action[®] Heel & Elbow Pads

Action's Heel & Elbow pads are recommended for people "at risk" or diagnosed with compromised skin disorders. Pad use allows clients taking medication that affects skin integrity, or who are prone to skin tears or abrasions, to continue daily activities or sleep comfortably.

Item#	Size
20401	.95 cm thick Heel & Elbow Pad (sold as pr)
20401A	.48 cm thick Heel & Elbow Pad (sold as pr)



Action[®] Donut Head Pads

Model 40201

20.5 OD x 7.5 ID x 4.5 cm

Designed to protect and cradle the patient's head.




Action[®] Horseshoe Head Pads

Model 40204X

20.5 OD x 7.5 ID x 2.8 cm

A low profile head pad to protect the patient's ear when lying in the lateral position.

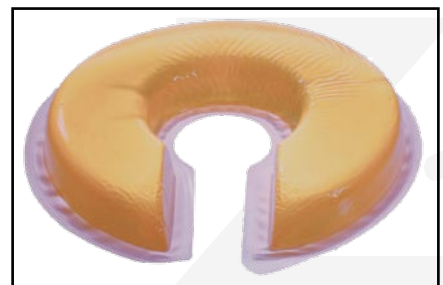


Action[®] Transfer Bench Pad

Model TB241678

61 x 41 x 2.2 cm

To reduce pressure when using a transfer bench for extended period. Custom sizes available.

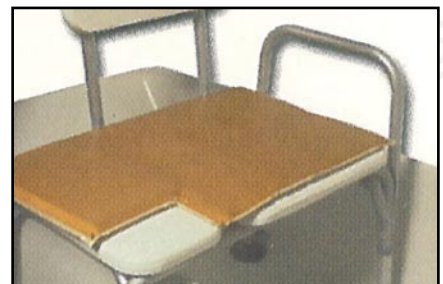

The Action® Professional™ Cushion



The Professional provides exceptional pressure relief for users at risk of skin breakdown or with a history of pressure ulcers. A full 4.1 cm thickness of Akton® polymer provides effective weight distribution on a contoured or flat surface. The Professional provides an additional 50% cushioning against pressure, shear, and vibration. This is significant for those active clients using a power wheelchair who have increased concerns about bottoming-out and the potential effects of increased temperatures (due to foam insulation properties). Custom sizes available upon request.

Professional Cushion with Basic Cover

4.13 cm AKTON polymer

Item #	Size, cm	Weight
5100-2	40.6 x 40.6	6.8 k
5105-2	45.7 x 40.6	7.2 k
5106-2	40.6 x 45.7	7.2 k
5108-2	45.7 x 45.7	9.0 k
5118202	45.7 x 50.8	9.9 k
5120182	50.8 x 45.7	9.9 k
512020	50.8 x 50.8	9.9 k

Standard Size Replacement Covers

COV51 Basic Cover

COI51 Incontinent Cover

The Action® Pilot™ Cushion



The Pilot is a low-profile cushion that provides exceptional pressure relief for most users with a pressure ulcer history. 2.54 cm Akton® polymer provides immersion for bony prominences and protects clients with fragile skin integrity. This multi-use cushion is recommended for use on manual and power wheelchairs, truck or taxi seats, or for other uses requiring significantly reduced high-pressure areas, vibration, and shear.

The low profile allows users to effectively reach the floor for more efficient foot-assisted propulsion. It also conforms to pre-contoured foam or rigid cushions to maximize stability and positioning when required.

Custom sizes available upon request.

The Basic Cover is sold with the Pilot. Incontinent and replacement covers are also available.

Pilot Cushion with Basic Cover

2.54 cm AKTON polymer

Item #	Size cm	Weight
9000-2	40.6 x 40.6	4.5 k
9006-2	40.6 x 45.7	4.9 k
9005-2	45.7 x 40.6	4.9 k
9008-2	45.7 x 45.7	5.4 k
9018202	45.7 x 50.8	6.3 k
902016	50.8 x 40.6	5.4 k
9020182	50.8 x 45.7	6.3 k
902020	50.8 x 50.8	6.8 k

Standard Size Replacement Covers

COV90 Basic Cover COI90 Incontinent Cover

The Action[®] Commuter[™] Cushion

The Commuter offers basic pressure and shear reduction for low risk users at an affordable price. Combined 1.27 cm Akton[®] polymer and 4.45 cm foam are sealed in a waterproof film for incontinence protection.

Commuter Cushion with Basic Cover

1.27 cm AKTON polymer with 4.45 cm foam

Item #	Size, cm	Weight
8500-2	40.6 x 40.6	2.2 k
851618	40.6 x 45.7	2.7 k
8505-2	45.7 x 40.6	2.7 k
8508-2	45.7 x 45.7	2.7 k
851820	45.7 x 50.8	3.1 k
8520182	50.8 x 45.7	3.1 k
852020	50.8 x 50.8	3.6 k

Standard Size Replacement Covers

COV85 Basic Cover



The Action[®] Centurian[™] Cushion

The Centurian offers pressure and shear reduction for low-to-moderate risk users. Combined layers of 1.6 cm Akton[®] polymer and 1.6 cm foam are sealed in waterproof film for incontinence protection.

Offering significant comfort for active clients who need a light to medium weight cushion with low maintenance, the Centurian is also the first choice of long term care facilities for decubitus prevention and treatment. Custom sizes available upon request. Maximum size of 60.9 x 60.9 cm in 2.54 cm increments.

Centurian Cushion with Basic Cover

1.6 cm AKTON polymer with 1.6 cm foam

Item #	Size cm	Weight
5200-2	40.6 x 40.6	2.7 k
5206-2	40.6 x 45.7	3.1 k
5205-2	45.7 x 40.6	3.1 k
5208-2	45.7 x 45.7	3.6 k
5218202	45.7 x 50.8	3.6 k
5220182	50.8 x 45.7	3.6 k
522020	50.8 x 50.8	3.6 k

Standard Size Replacement Covers

COV52 Basic Cover

COI52 Incontinent Cover



The Xact® Classic™ Cushion



Standard Size Replacement Covers

COBX BODY SMART

CORX Incontinent Cover

The Xact Classic offers outstanding pressure distribution and anti-shear properties. The 1.4cm Action® cube pad and laminated foam base help the user maintain a comfortable, therapeutic position. Rehab therapists most often recommend the Classic for users requiring forward and lateral positional support. The recent re-design that incorporates viscoelastic pink foam offers improved support for gluteal areas and bony protuberances. Custom sizes available in 2.54 cm increments with standard pommel and side rails. Maxium size is 60.9 x 60.9 cm. A Body Smart® Cover is sold with the Classic. Incontinent and replacement covers are also available.

Classic Cushion with BODY SMART® Cover

Item #	Size, cm	Side Rails	Weight
X1010	25.4 x 25.4	2.54 cm	1.3 k
X1012	25.4 x 30.5	2.54 cm	1.3 k
X1210	30.5 x 25.4	2.54 cm	1.3 k
X1212	30.5 x 30.5	2.54 cm	1.8 k
X1214	30.5 x 35.6	2.54 cm	1.8 k
X1412	35.6 x 30.5	2.54 cm	2.2 k
X1414	35.6 x 35.6	2.54 cm	2.7 k
X1416	35.6 x 40.6	2.54 cm	3.1 k
X1614	40.6 x 35.6	2.54 cm	3.1 k
X1616	40.6 x 40.6	3.81 cm	3.1 k
X1618	40.6 x 45.7	3.81 cm	4.5 k
X1620	40.6 x 50.8	3.81 cm	4.5 k
X1717	43.2 x 43.2	3.81 cm	4.9 k
X1816	45.7 x 40.6	3.81 cm	4.5 k
X1818	45.7 x 45.7	3.81 cm	4.9 k
X1820	45.7 x 50.8	5.08 cm	4.9 k
X2016	50.8 x 40.6	5.08 cm	4.5 k
X2018	50.8 x 45.7	5.08 cm	4.9 k
X2020	50.8 x 50.8	5.08 cm	5.8 k

Cushion Cover Information

Basic Cover - The Basic Cover is made of a durable stretch knit material which is essential to utilize the protective cushion properties. Moist skin can significantly compromise skin integrity, contributing to decubitus ulcers. This material allows voids to pass through and optimize air exchange, keeping clients dryer as it wicks moisture away.

This is a vital issue for the professional client with limited time to address personal needs, and for long term care residents. This cover is standard on the Commuter™, Centurian™, Pilot™ and Professional™ cushions.

Body Smart® Cover - The BODY SMART Cover is designed to assist in the proper alignment and positioning of the client. With its breathable double layer, the cover provides support in the ischial area against the effects of shearing, while the front of the cover offers friction to maintain hip position and encourage stability. This cover is standard on the Xact® Line.

Incontinent Cover - This Low Shear Cover is designed to protect the user from full contact with voids by allowing them to flow through to the moisture-resistant layer. This style of Incontinent Cover is intended for the Centurian™, Pilot™, and Professional™.

Reverse Dartex® Cover - Designed for the XACT Cushion series, this cover offers excellent moisture-resistant options for incontinent use. This style resists voids and is easy to clean. The Reverse DARTEX cover has a very smooth, non-fabric finish.

Shear Smart® Cover - The SHEAR SMART Cover provides significant relief from the effects of shearing and is an effective moisture barrier against incontinence. It offers a complete double layer of low shear material to allow for maximum movement on the pad. When combined with the Cube Pad, it significantly decreases the consequences of shearing while providing comfort from the effects of pressure.

The XACT[®] Lite

Action Products XACT Cushion range continues to expand with the addition of the XACT LITE, a cushion for the active and mobile patient, including those who require sports-friendly cushions, and can self-propel, transfer and shift.

The new design of the AKTON[®] polymer reduces shear and pressure and provides maximum surface flexibility with no hammocking effect. Additional features include:

- Light leg troughing which improves stability and posture without interfering with transfers
- Low profile design which helps foot propulsion
- Pre-ischial ridge to help maintain pelvic neutrality
- Dual-density foam base with firm rear corners to improve stability, and soft front corners to aid transfers
- Cushion weight of 1.45kg, and height of 5.4 cm

The XACT LITE has been tested for a weight capacity of 114kg.

SIZE RANGE

XL 1616	40.6 X 40.6 cm	XL 1820	45.7 X 50.8 cm
XL 1618	40.6 X 45.7 cm	XL 2016	50.8 X 40.6 cm
XL 1620	40.6 X 50.8 cm	XL 2018	50.8 X 45.7 cm
XL 1816	45.7 X 40.6 cm	XL 2020	50.8 X 50.8 cm
XL 1818	45.7 X 45.7 cm		

The XACT[®] Soft

The XACT SOFT provides excellent comfort to those who have partial or full sensation, and are mid to high risk ulcer candidates.


The unique construction of the XACT SOFT is the result of Action's on-going research in pressure care, and provides:

- Postural stability
- A sacral dish conforming to the user
- An Action[®] CUBE insert minimising shear and cushioning the IT area
- A Body Smart[®] cover with a gripping front layer to help prevent slouching
- Easy handling due to its light weight of 1.7kg, Height is 11cm.

The XACT SOFT has been tested for weight capacity of 160kg.

SIZE RANGE

XS 1616	40.6 X 40.6 cm	XS 1820	45.7 X 50.8 cm
XS 1618	40.6 X 45.7 cm	XS 2016	50.8 X 40.6 cm
XS 1620	40.6 X 50.8 cm	XS 2018	50.8 X 45.7 cm
XS 1816	45.7 X 40.6 cm	XS 2020	50.8 X 50.8 cm
XS 1818	45.7 X 45.7 cm	(CUSTOM SIZES NOT AVAILABLE)	

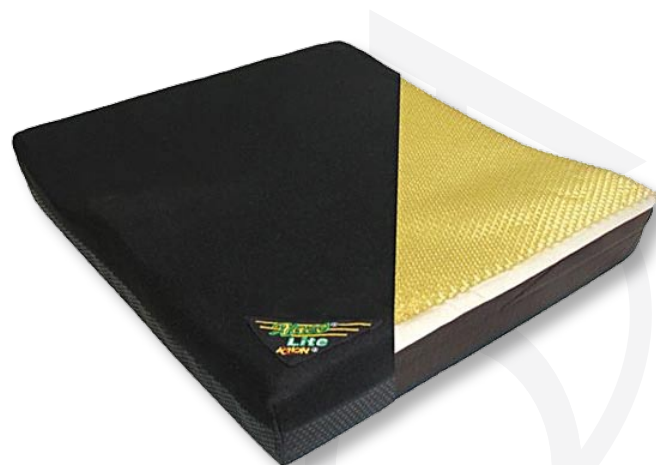


Xact[®] Individual Cushion

The Individual allows a seating specialist or experienced therapist to custom-build an individual seating cushion tailored to wheelchair user's specific seating needs. By using the various XACT Positioning Assessment Kit component pieces, a therapist can build a personalized seating product that provides comfort, low shear, corrective positioning, and high-quality pressure management.

The therapist can actually try various seating intervention strategies before deciding on a final configuration. Individual cushions are customized and ordered by the therapist and assembled into a final product by Action Products, Inc. Thousands of unique seating combinations are available.


The Body Smart[®] Cover is sold with the Individual.

Incontinent and replacement covers are also available.

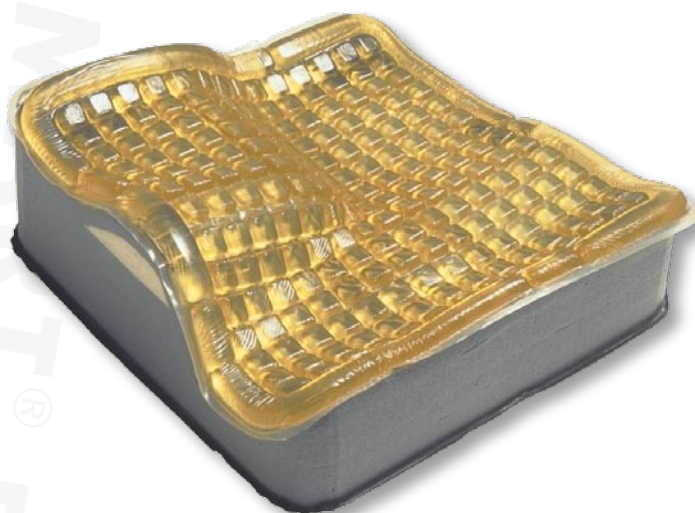

Individual Components

- A.** Action[®] Cube Pad
- B.** Adductor Rails
- C.** Left Obliquity
- D.** Left Contracture
- E.** Pommel
- F.** Anti-Thrust
- G.** Individual Cushion Base
- H.** Right Contracture
- I.** Right Obliquity

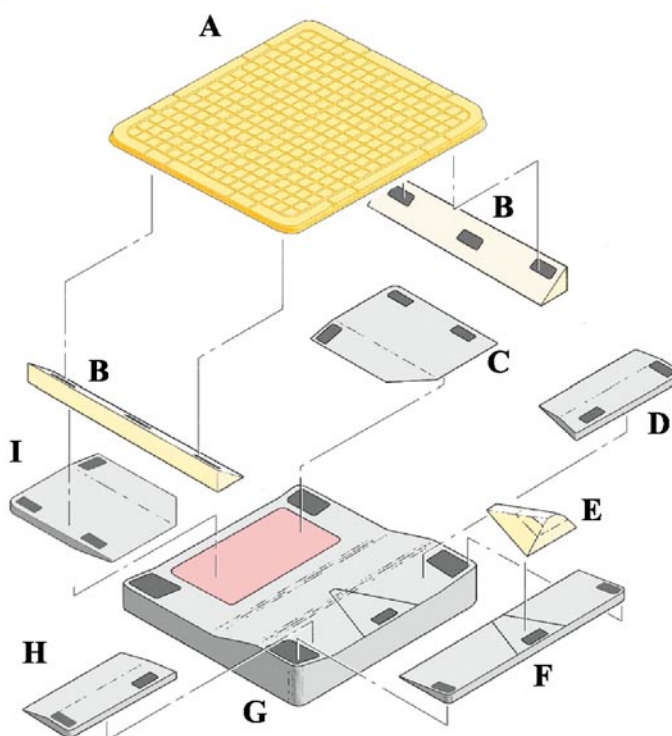

For detailed information request separate brochure.

Shear Only
Abduction
Adduction
Atrophy
Existing Ischial Injury
Existing Sacral Injury
Foot Propel
Impaired Muscle Control
Skeletal Deformities
Pelvic Stability
Positioning
Risk-Low
Risk-Medium
Risk-High
Windswept
Tilt/Wheelchairs
Active Users
Sports

Standard Sizes in CM


Commuter™

Custom Sizes Available

CenturianTM

Custom Sizes Available

Professional™

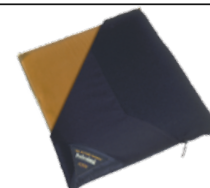

Custom Sizes Available

Pilot™



Custom Sizes Available

Xact®



Custom Sizes See Xact Individual

Mattress Overlay



Mattre



Heel & Elbow Pads



Transfer Bench Pad





Standard Size Replacement Covers

#COG SHEAR SMART

(Specify size when placing order)

The Shear Smart® Pad

The SHEAR SMART pad, combined with the breathable double-layer SHEAR SMART cover, significantly reduces the effects of shear while providing superb user comfort. The cube construction of the pad surface allows air circulation, provides cushioning from mild impact, heat dissipation, and the comfort of sitting on 1.4 cm Akton® polymer.

This cushion will conform to contoured bases for additional support when required. Custom sizes available upon request. Maximum of 60.9 x 60.9 cm in 2.54 cm increments.

The SHEAR SMART Cover is sold with the SHEAR SMART Pad. Replacement covers are available.

SHEAR SMART Pad with SHEAR SMART Cover 1.4 cm Action® cube pad with cover

Item #	Size, cm	Weight
CG1616	40.6 x 40.6	1.81 k
CG1618	40.6 x 45.7	2.27 k
CG1816	45.7 x 40.6	2.27 k
CG1818	45.7 x 45.7	2.72 k
CG1820	45.7 x 50.8	2.72 k
CG2018	50.8 x 45.7	2.72 k
CG2020	50.8 x 50.8	3.17 k

Twister® Wheel Chair Back

The Twister Back is a hybrid back support using hook and loop technology to support users with normal to severe kypho-scoliotic vertebral misalignment. Akton® polymer cubes horizontally line the back of the cover and are designed to accommodate deformity while ensuring comfort and pressure relief. The cube construction make the product exceptionally lightweight and allows air to flow between rows, keeping it cool. This lightweight design makes this product ideal for active users.

Available in two styles, Short and Tall, the backs are for both cane back and rigid frame chairs.

When ordering, please specify if product is for chair with or without push handles.

TWISTER Tall Back: for chair backs 43cm or higher

Item#	Fits these dimensions
TCNT	38-41cm Width x 43-53 cm Height
TCAT	43-46cm Width x 43-53 cm Height
TCMT	48-51cm Width x 43-53 cm Height
TCWT	53-58cm Width x 43-53 cm Height

Specify PHT for Cane with push handles or WPHT for Cane without push handles.

TWISTER Short Back: for chair backs 41cm or lower

Item#	Fits these dimensions
TCNS	38-41cm Width x 28-41 cm Height
TCAS	43-46cm Width x 28-41 cm Height
TCMS	48-51cm Width x 28-41 cm Height
TCWS	53-58cm Width x 28-41 cm Height

Specify PHS for Cane with push handles or WPHS for Cane without push handles.


The Action® Adaptive Cube Pad

Action's Adaptive Cube Pads also have multifunctional applications for traditional and non-traditional seating challenges. Pads are designed to be cut and used for spot protection in products such as arm pads, foot plates, head/leg rests and splints/braces. Pads available in 1.4 cm thickness. The Adaptive Cube Pad is not sold with a cover. If a cover is desired, please order a Shear Smart® Pad as described . Custom sizes available upon request.

Maximum of 60.9 x 60.9 cm in 2.54 cm increments.

Item #	Size	Weight
CU1616	40.6 x 40.6 x 1.4 cm	2.27 k
CU1618	40.6 x 45.7 x 1.4 cm	2.27 k
CU1816	45.7 x 40.6 x 1.4 cm	2.27 k
CU1818	45.7 x 45.7 x 1.4 cm	2.72 k
CU1820	45.7 x 50.8 x 1.4 cm	2.72 k
CU2018	50.8 x 45.7 x 1.4 cm	2.72 k
CU2020	50.8 x 50.8 x 1.4 cm	3.17 k

The Adaptive Cube Pad

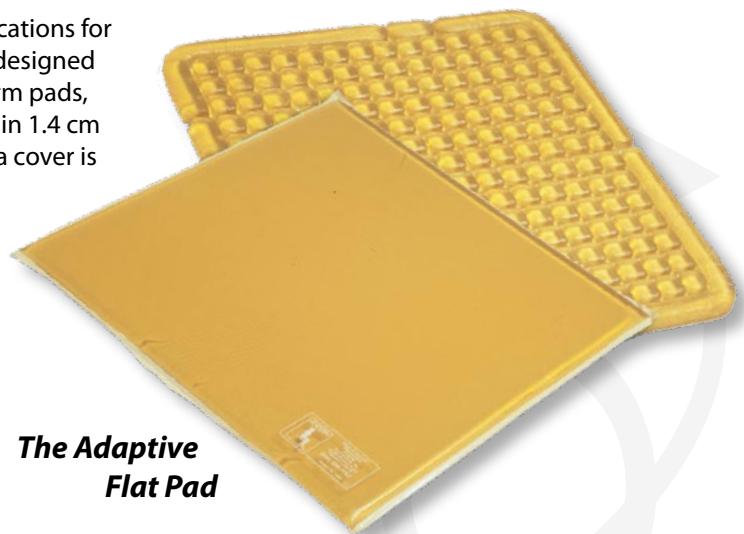

The Adaptive Flat Pad

The Action® Adaptive Flat Pad

Action's Adaptive Flat Pads have multifunctional applications for traditional and non-traditional seating challenges. Pads can be cut and used for spot protection on products such as arm pads, foot plates, head/leg rest and splints/braces. Pads are available in .635 cm, 1.3 cm, or 1.59 cm thickness. The Adaptive Flat Pad when used as a cushion is intended for comfort only, not medical intervention. Note: If 1.3 cm or 1.59 cm pad is used for a seating cushion, pad must be used with an Action® brand cover for warranty to be in force.

The Adaptive Flat Pad comes standard without a cover, however the **Basic Cover** is available for all models except the 561818.

Item #	Size, cm	Weight
561818	45.7 x 45.7 x .635 cm	1.81 k
5302	40.6 x 40.6 x 1.3 cm	2.27 k
5303	45.7 x 40.6 x 1.3 cm	2.72 k
541818	45.7 x 45.7 x 1.3 cm	2.72 k
6500	50.8 x 50.8 x 1.3 cm	3.20 k
5300	40.6 x 40.6 x 1.59 cm	2.72 k
5301	45.7 x 40.6 x 1.59 cm	3.17 k
531818	45.7 x 45.7 x 1.59 cm	3.62 k

Optional Covers:

Covers are only available for 1.3 cm and 1.59 cm Adaptive Flat Pads. To order, use Item #COV53 followed by the pad dimensions.

Eg: COV531616 will fit a 40.6 x 40.6 cm Adaptive Flat Pad.

Hand Exerciser

Made of 100% Akton® polymer, these flexible cylindrical exercisers provide uniform resistance for individual finger or gross grasping exercises.

The polymer provides unlimited degrees of resistance directly proportionate to the force placed on it. The Hand Exerciser is perfect for increasing hand strength and coordination of people of all ages.

Repair Kits

Use either the Patch or the Roll for quick and easy resurfacing of small cuts or punctures on any Action® pad. Each Repair Kit contains full instructions for use.

Item #	Description
KIT002	Resurfacing Patch 10x10cm
KIT003	Resurfacing Roll 91x10cm



Item #	Description
70001	Soft Hand Exerciser
70003	Firm Hand Exerciser

ACTION® ACTION®, AKTON®.

BODY SMART®, CENTURIAN™, PILOT™, PROFESSIONAL™, COMMUTER™,
XACT®, SHEAR SMART®, TWISTER® ARE EITHER REGISTERED TRADEMARKS
OR TRADE MARKS OWNED BY ACTION PRODUCTS, INC. U.S.A.

DARTEX IS A REGISTERED TRADEMARK OF DARTEX COATINGS, INC. U.S.A.

Your Distributor is:



SURGICAL AND MEDICAL SUPPLIES PTY LTD

16 Kensington Road Rose Park SA 5067

Telephone: (08) 8332 1666 Fax: (08) 8332 5132

Email: sales@surmed.com.au

Website: www.surmed.com.au

Imported and Distributed in Australia by:



EDWARDS
M E D I C A L

15 Bellingham Street, (P.O. Box 3085)
Narellan NSW 2567

P. (02) 4647 1444

F. (02) 4647 9871

Email: info@edwardsco.com.au

© 2007