

Stretcher Mattress

High grade mattress with excellent durability and pressure relieving comfort for patient stretchers.

Features

- Thicknesses that aid egress and imaging
- Premia mattress - 3 effective, layered foams
- Protran mattress - single foam layer
- Welded 'double coated' stretch cover
- Non-slip base design
- Permanent care label
- New technology fluid proof zip

Product Code

- 100mm Premia Mattress; PT-CM-100-PRP
- 150mm Premia Barituff Mattress; PT-CM-150-PRB
- 100mm Protran Mattress; PT-CM-100-PTP
- 75mm Protran Mattress; PT-CM-75-PTP

Mattress Size





- Standard; 680 W x 1900mm L
- Narrow; 570mm W x 1900mm L
- Wide; 830mm W x 1900mm L
- Barituff; 980mm W x 1900mm L

Mattress Shape

- Square Ends (electric stretcher)
- Truncated End (hydraulic stretcher)

Summary

The **Contour Premia** mattress has been designed for an increased patient weight capacity and comfort by using three quality foam layers, a high density layer to assist in preventing patients from bottoming and a slow recovery, memory foam upper layer to reduce pressure. It is supplied with a quality welded, stretch and fluid proof cover with permanent label to assist in caring for the mattress. The Contour Protran mattress option is only recommended for short stay patients with a single medium density foam section to facilitate special patient transfers or for clear imaging.

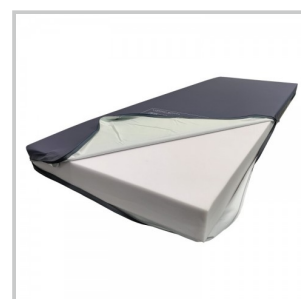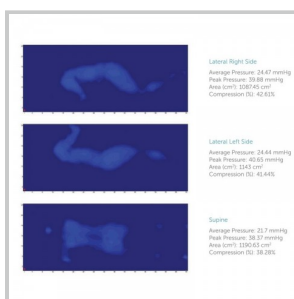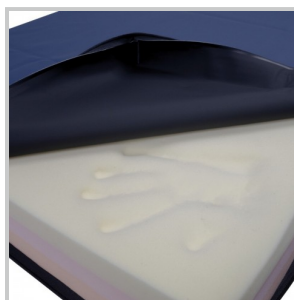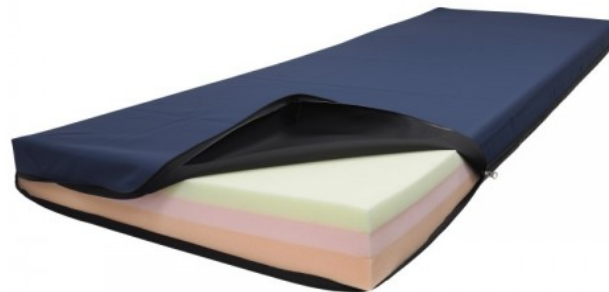

12/11/2019