

Reynard's Premier detergent & disinfectant wipes.



Surgical & Medical Supplies

Stay one wipe ahead of germs.

Clean and disinfect equipment and surfaces simply and effectively whilst saving time and money.

Reynard's Premier detergent and disinfectant wipe is a proven and single step solution to cleaning and disinfecting surfaces equipment. Our new gentle formula is optimised for use in both clinical and non-clinical settings on surfaces and equipment whilst proven to kill at least 99.999% of pathogens and is tested and ARTG listed for use against SARS CoV 1&2. AUSTL337692.



Reynard is reassurance.



Reynard
Health Supplies

Reynard's Premier detergent & disinfectant wipes.


The cloth that is a cut above.

Our Premier wipes are gentler on materials and have a newly developed soft apertured polyester substrate that boasts superior formulation retention, and extremely good traction to pick up dirt, bacteria and other pathogens. Most other surface disinfecting products are made from thinner, thermal bond polyester substrates that push the dirt or bacteria around, rather than pick it up.

This provides significantly higher levels of performance especially if you are using a 1 step cleaning and disinfecting protocol. Combined with a softer formulation that will not dry or damage surfaces and is safe to use on surfaces and equipment, Reynard Premier disinfectant wipes provide a safe, easy to use and cost-effective infection control solution.

The All-in-1 wipe that has best in class performance for cleaning and disinfection.

Reynard Premier Disinfectant Wipe

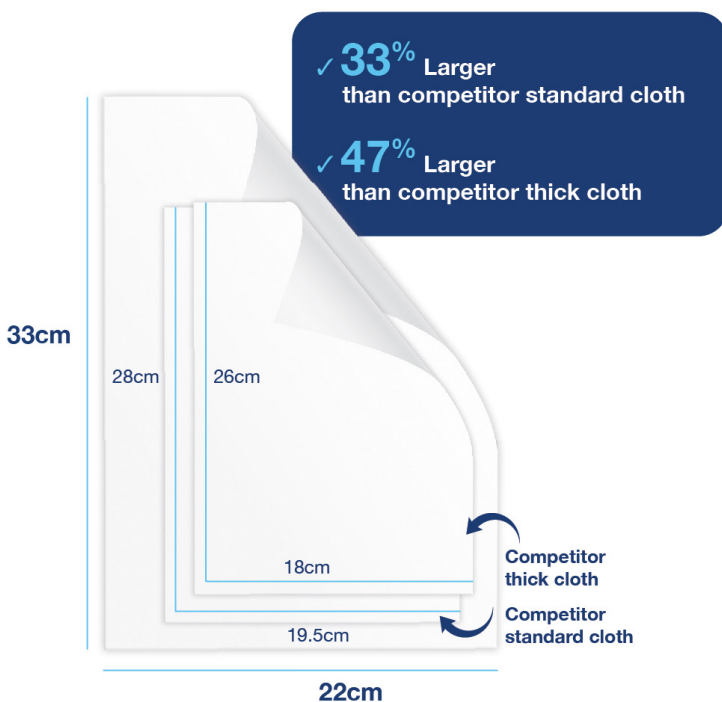

Our winning formula.


Biodegradable non-ionic surfactant with excellent non-streaking and superior wetting properties to provide the wipe with an effective cleaning action and enable the liquid to effectively wet out surfaces.

It incorporates a Fungicide and Chelating agent to enhance the biocidal activity of the wipe. Our proprietary active disinfecting formula that contains three active biocidal agents for maximum effectiveness and adaptability.

While germs aren't safe, your equipment is.

Reynard's Premier wipes were developed with key clinical equipment providers to ensure an unparalleled range of device and material compatibility. The wipes will not dry or damage surfaces, and are safe to use on plastics, polymers, aluminium, stainless steel, powder-coated surfaces, vinyl and leather. Tested and safe to use on and around equipment. The formula was developed so that the wipes would clean hard surfaces & equipment without causing harm.

✓ Reynard Premier Disinfectant Wipe is a **28%** thicker cloth

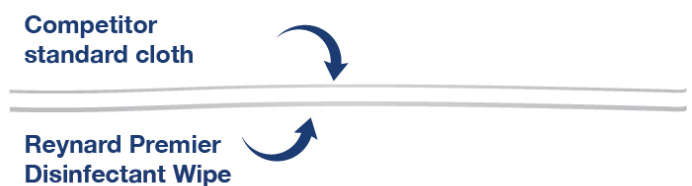



Kills 99.999% of pathogens after just 60 seconds, bacteria, fungi, viruses including SARS-CoV-1&2 and MERS-CoV.

Reynard Health Supplies.

Reynard Health Supplies was established in New Zealand in 1999 and today is the number one supplier of surface cleaning & disinfection wipes, specialist patient care wipes and IV Therapy products to hospitals in New Zealand and Australia.

Reynard Health Supplies is owned by a dedicated group of Australian and New Zealand shareholders who are all actively and directly involved in the day to day operation and development of the business worldwide.

Reynard is reassurance.

At Reynard, we strive to provide products and support services that provide you with the reassurance you need to be able to focus on optimum patient outcomes. This reassurance stems from:

- Quality & Range - a wide range of quality products to suit every need;
- Cost savings - Cost effective solutions;
- Reliability - A reliable, proven supply chain;
- Support - Technical, educational and compliance support services;
- Research & Development - Ongoing product development; and
- Flexibility - The willingness to design new products and supply chain options.

Reynard is reassurance.



Product technical summary and ordering information.

RHS216 Premier disinfectant wipe

- 50 wipes per pack
- 12 Packs Per Carton
- 33 x 22cm (25 to 50% larger than most competitors)
- Thicker, soft apertured polyester substrate
- Registered Disinfectant with claims

RHS218 Premier disinfectant Wipe

- 100 Wipes Per Pack
- 6 Packs Per Carton

Related products.

- RHS201, 203 and 206 Neutral Detergent Wipes
- RHS202, 204 and 205 Surface Disinfectant Wipes
- RHS402 and 409 IPA Surface Wipes

Testing Data.

All claims in accordance with TGO 104 Order 2019

Effective against a wide range of pathogenic micro-organisms including bacteria, viruses such as SARS-CoV-1&2 and MERS-CoV virus and fungi.

PATHOGEN	EN Test	Time
Bacteria		
Kills at least 99.999% of pathogens. Effective from 60 seconds		
Enterococcus Hiera	16615	1 min
Staphylococcus Aureus	16615	5 min
Fungi		
Aspergillus Niger	13624	5 min
Candida Albicans	13624	1 min
Candida Albicans	16615	5 min
Candida Auris	13624	1 min
Viruses		
Effective against enveloped viruses		
SARS-CoV-1&2 and MERS	14416	1 min

Also tested to European standards EN1276, EN13727, EN14348, EN13704, EN1650, EN14675, EN14476, EN14651, EN1656, EN16615



SA Distributor:

Australia

44 Hampden Road
 Artarmon, NSW, 2064
 enquiries@reynardhealth.com.au
 1300 667 199

New Zealand

PO Box 8470
 Havelock North
 Sales@reynardhealth.com
 0508 494 737

International

PO Box 8470
 Havelock North
 Sales@reynardhealth.com
 +64 (0)6 650 0709



SURGICAL AND MEDICAL SUPPLIES PTY LTD
 16 Kensington Road Rose Park SA 5067
 Telephone: (08) 8332 1666 Fax: (08) 8332 5132
 Email: sales@surmed.com.au
 Website: www.surmed.com.au